



BioTemps® NOW

“The Temporary You Admire”



Affordable PMMA material designed for **esthetic, in-office-milled temporaries**

- Recommended for temporary full-coverage anterior and posterior crowns
- Mills in under 10 minutes*
- Available exclusively for Glidewell milling systems
- Low water sorption means long-lasting strength and shade retention



PATIENT-PLEASING RESULTS



NO OVEN NEEDED



RETAINS GINGIVAL POSITION

Before



After



This patient presented with failing restorations and recurrent decay on teeth #12–15. Because of the number of units involved, BioTemps® NOW provisional restorations were milled and seated for a durable and comfortable solution while awaiting final treatment.



For our cosmetically-driven cases, BioTemps NOW is an excellent option to meet the demand for esthetics while the final restorations are being fabricated. Furthermore, the pre-treatment scanning with the fastdesign.io software allows the temporaries and final restorations to mimic our diagnostic models and wax-ups.

- Joshua Prentice, DDS | Bend, OR

For more information
glidewelldental.com/milling-blocks ■ 844-949-7184



BENEFITS

- Give your patients the ability to preview their final results
- Deliver excellent abrasion resistance for up to six months of wear time



Sold in 5-block packs

READY TO DELIVER IN UNDER 15 MINUTES



MILL

+



FINISH/POLISH

+



NO OVEN
REQUIRED

=



TOTAL TIME

SPECIFICATIONS

MATERIAL

- PMMA
- Flexural strength: ≥ 50 MPa

IDEAL FOR


- Temporary full-coverage anterior and posterior crowns

SHADES


- 6 VITA® Classical Shades, 1 Bleach Shade: A1, A2, A3, A3.5, B1, C2, BL

EASY PREPARATION

POSTERIOR

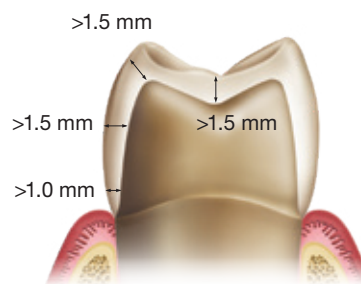

ANTERIOR

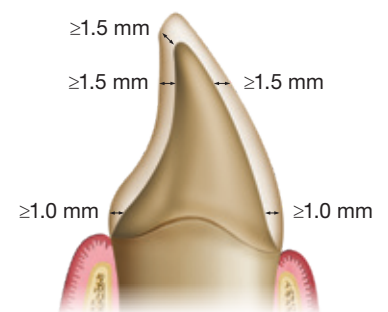

- Wall thickness should be at least 1.5 mm (1.0 mm at margins)
- Margins should be prepared with a deep chamfer or rounded shoulder

*Milling times vary based on prescription and mill used. Average milling time is shown. VITA is a registered trademark of VITA Zahnfabrik.

MKT-10465_4 PK-281915-022720